


Let's Sing

Shower Stool

Article no. LSG101 White

Article no. LSG201 Black

Let's Sing is a new shower stool in a modern design intended for use in the shower or the bathroom. The stool features a beautiful and stable construction for safe and secure use. The legs are easy to adjust to desired height with one easy motion. A wire basket is available as an accessory that should be placed under the seat for storing shower articles or other products.

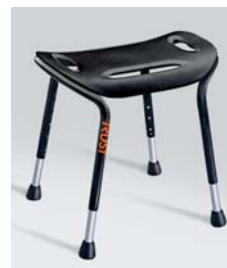

Let's Sing is adjustable from 45 to 56 cm. The height is adjusted by pressing the buttons on the legs and pulling the lower legs to the desired positions, after which the buttons will be locked in place in the adjustment holes.

Important: Always adjust the legs to the same height to avoid falling accidents.

Technical data:

Material:
 Seat ABS
 Tube fram aluminum
 Feet TPE
 Height adjustable from 45 to 56 cm
 Weight 2,4 kg
 Tested according to ISO EN 12182:2012
 Max user weight 150 kg

Warranty:
 Two years warranty

